



Barrier Free Architecturals
480 E Franklin Rd, #108, Meridian, ID 83646
1(877) 717-7027 or visit: www.barrierfree.org
EMAIL: info@barrierfree.org

Installation Instructions

Semi-permanent Threshold Adaptors: TA-S, TA-SW, TA-SM

TAS

Supplies provided:

Threshold Adaptor
Silicone (100%) Bath tub Caulk, white
Instructions

Tools and supplies needed:

Pencil
Masking tape (optional)

To install the Adaptor:

1. Put the Adaptor in place and trace a pencil line around it on the shower bottom. See Figure 1
2. Place a 1/4" bead of silicone caulk just inside the pencil line.
3. Carefully put the threshold in position. See Figure 2
4. Smooth out the excess silicone with a wet finger. Masking tape positioned 1/4" from the pencil line on the shower and 1/4" from the edge of the Adaptor can help here. See Figures 2 and 3
5. Remove the masking tape and touch up with a wet finger to feather in the silicone.
6. Do not disturb for 12 hours.

To remove the Adaptor:

1. Carefully cut through all the caulk with an Exacto-type knife, holding the blade flat.
2. Pivot the threshold into the shower and pull it up.
3. Carefully remove the caulk with a single edged razor blade and "Goof Off" Remover (mfd. by Guardsman, available at most hardware stores)

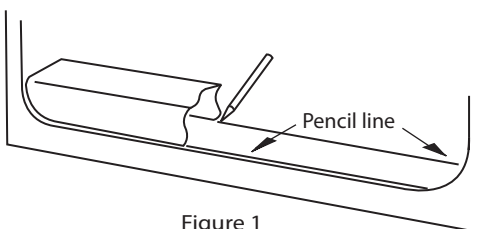


Figure 1

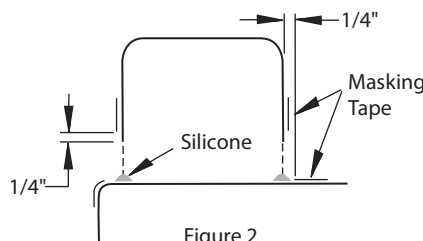


Figure 2

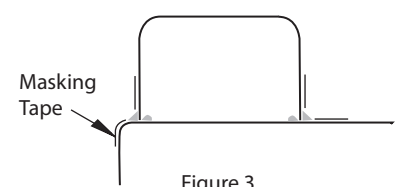


Figure 3